




HNPS Advances in Nuclear Physics

Vol 11 (2002)

HNPS2000 and HNPS2002

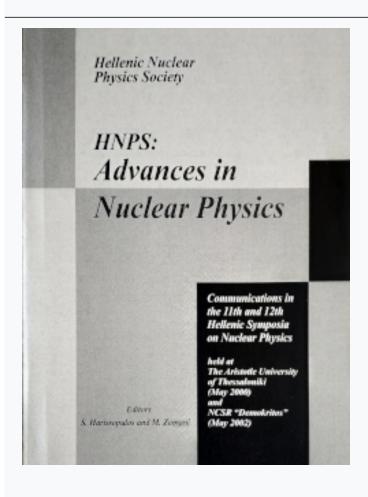

Analysis of Ancient Glass Using Ion Beams and Related Techniques

S. Kossionides, M. Kokkoris, A. G. Karydas, T. Paradellis, G. Kordas, G. Moraitou

doi: 10.12681/hnps.2226

To cite this article:

Kossionides, S., Kokkoris, M., Karydas, A. G., Paradellis, T., Kordas, G., & Moraitou, G. (2019). Analysis of Ancient Glass Using Ion Beams and Related Techniques. *HNPS Advances in Nuclear Physics*, *11*. https://doi.org/10.12681/hnps.2226

Analysis of Ancient Glass Using Ion Beams and Related Techniques.

S. Kossionides¹, M. Kokkoris¹, A. G. Karydas¹ and T. Paradellis¹ G. Kordas², G. Moraitou^{2a}

¹Material Analysis Laboratory, Institute of Nuclear Physics, NCSR Demokritos ²Sol-Gel Laboratory, Institute of Materials, NCSR Demokritos ^aMinistry of Culture

Abstract

Ion Beam Analysis (IBA) techniques are established as an indispensable tool for material analysis. They are very powerful when applied to the examination of uniform layers. In recent years they have been also used for the quality control of implanted layers, which improve the mechanical or chemical properties of common materials. The use of IBA techniques for the study of corrosion is impeded by the non-uniformity of such layers, which renders the analysis difficult. In the Material Analysis Laboratory of the Institute of Nuclear Physics, a combination of experimental techniques and analytical algorithms has been developed, which allows the non-destructive examination of corrosion layers on cultural artifacts.

Ancient glass fragments from Kenchreai, Greece, were analyzed with ion beam techniques, in order to examine the applicability of the latter to the characterization of the surface. The relative distribution of heavy components was determined by X-Ray Fluorescence (XRF). Nuclear reactions were used to determine the depth distribution of H, O, C and Si. The Hydrogen distribution rises sharply from the surface, reaching a maximum of 9×10^{22} at/cm³ at 0.5 μ m and then decays smoothly to a depth of approximately 2 μ m. Examination of the same samples with α -particles, indicates a depletion of Silicon in the same surface region, while the Oxygen distribution is uniform. A detailed analysis, using the 28 Si(p, p) reaction, which resonates at 2.095 MeV, confirms the depletion and other non-uniformities in the Silicon distribution.